

# BROEN-LAB EMERGENCY SHOWER SYSTEMS



reddot design award  
winner 2013



EMERGENCY  
SHOWER SYSTEMS

# BROEN-LAB

## TABLE OF CONTENT

Table of Content / Product Overview	2
Introduction / Customer Statements	3

### Emergency Showers & Eye Washes


Hand-Held Eye Showers & Accessories	4-5
Eye/Face Showers	6
Recessed Emergency Showers	7-8
Body Showers	9
Combination Showers	10
Accessories for Emergency Showers	11-12
Connections for Emergency Showers	13
Test Equipment for Emergency Showers / Eye Wash Bottles	14

### Heated Emergency Showers & Eye Washes

Introduction	15
Heated Combination Showers	16
Heated Body Showers	17
Heated Eye/Face Showers	18
Mobile & Portable Emergency Showers	18
Heated Tank Showers	19
Accessories for Heated Showers	20
General Information	21
International Standards	22-23
Technical Information	24
Notes	25-27

## PRODUCT OVERVIEW

- Please refer to the actual page for product versions

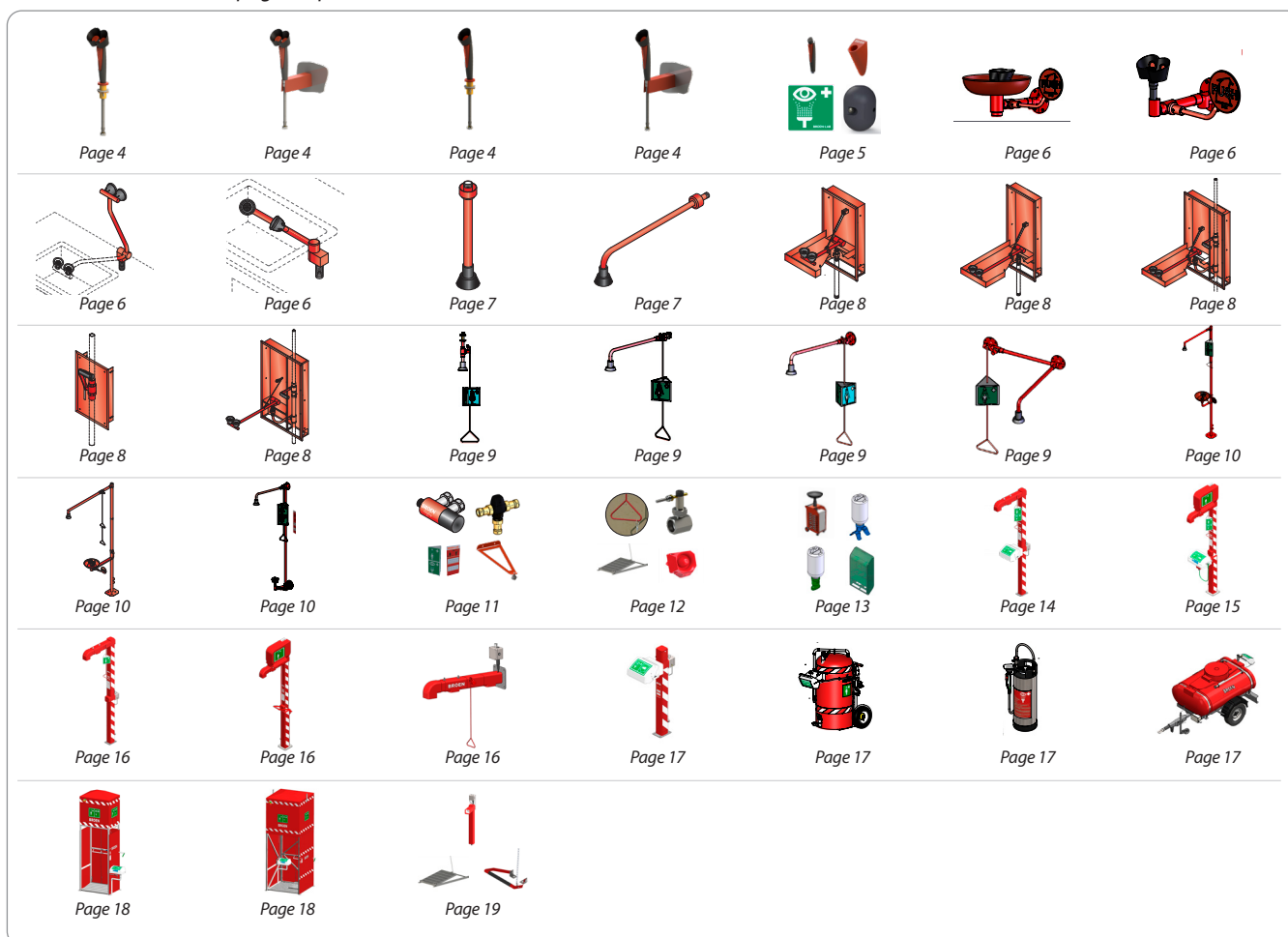

Our extensive knowledge allows us to continually develop and optimize our products so that they meet our customer needs and specifications. Safety, functionality and quality are an integral part of product design and flexibility.

## BROEN-LAB – A SAFE CHOICE

- Quality since 1948
- Complete product range
- Compliance with international standards
- Proven functionality world-wide
- Spare parts available for a guaranteed 10-year period



## CUSTOMER STATEMENTS



### Falck

#### IN THE SERVICE INDUSTRIES



"For more than four years Falck has had an excellent cooperation with BROEN-LAB. Our customers must be able to rely on Falck to guarantee the best solutions and consequently we only provide service for emergency showers that comply with the requirements of Danish Standard 15154 Parts I and II.

BROEN-LAB products are extremely easy to install and they are also clearly marked and easily serviced. Quite simply, they provide a better and safer solution for all parties."

Michael Hansen, Product Manager  
Falck



#### IN INDUSTRY



"At Carlsberg we do not accept a compromise with regard to quality. We have therefore replaced all our emergency shower systems by products from BROEN-LAB so that they comply with the EU norm EN 15154.

We want to look well after the people working for us and it is our distinct opinion that our workplace must be a place where you feel well and safe."

Jens Hansen,  
Coordinating Work Environment Representative  
Carlsberg Denmark



#### IN THE LABORATORY



"For quite some time our company has been in charge of the maintenance and service of installations at Akademiska Hus at Medicinberget in Gothenburg. This includes a number of laboratories which again means a considerable number of emergency shower systems supplied by BROEN-LAB.

We check the emergency showers regularly and are very content with BROEN-LAB's reliable design and construction which make them easy to inspect and test."

Anders Bergvall, Project Manager  
VVS-Teknik Mölndal



#### IN THE EDUCATIONAL SECTOR



"At DTU we teach more than 7,000 students within very different fields – from the metal industry to microbiology – and within each specific field unique safety aspects must be considered. Therefore it is essential for us that we have a cooperation partner who supplies a quality product and who can also advice and guide us on the specific risk profile of each classroom."

Bjarne Andersen,  
HVAC Manager at DTU Campus Service



#### IN THE MEDICAL INDUSTRY



"As an international advisor it is important for us to have a partner that complies with both national and international norms and standards – and who can document this by an independent third party certification. At the same time it is important to have a partner that can offer customized solutions supporting the architectonic expression of the project instead of colliding with it. Even though we are talking about emergency equipment, beautiful design is a plus."

Henrik Als Nielsen,  
Engineer, Ventilation & Building Utilities,  
NNE Pharma plan



## HAND-HELD EYE SHOWERS

BROEN-LAB has developed hand-held eye showers, which comply to EN15154-2 and ANSI Z358.1 norms and also offer many improved features.



### MATERIALS WITH MEDIA CONTACT

Plastics: PA, POM, TPE  
Metals: SS, AISI 316  
Rubber: EPDM

### SAFETY

- Improved trigger lock
- Improved activation of the shower
- BROEN-LAB recommends the installation of a Thermostatic Mixing Valve (TMV) with safety shower equipment for safe operation, see page 11



### OPERATION

- Stay-open dust cap
- Light weight
- Well-balanced construction
- Spray pattern



### RELIABILITY

- Easy-to-clean spray head
- Anti-limescale design and material
- Improved functionality
- Long lifetime



reddot design award  
winner 2013




### Eye Shower with Two 45° Heads



No. 27 10 009 10 05 011

Weight: Approx. 0.5 kg

Table mounted eye shower with two 45° angled heads, incl. 1.5m inlet hose. Trigger operated including a built-in check valve and FLOWFIX for regulation of water flow (14 l/min). Min. working pressure: 1,5 bar.

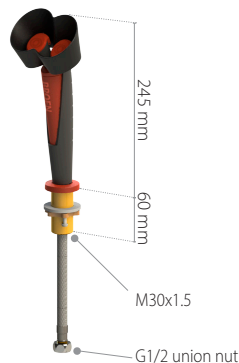


### Eye Shower with Two 45° Heads



No. 27 10 009 10 03 011

Weight: Approx. 0.5 kg

Wall mounted eye shower with two 45° angled heads, incl. 1.5m inlet hose. Trigger operated including a built-in check valve and FLOWFIX for regulation of water flow (14 l/min). Min. working pressure: 1,5 bar.

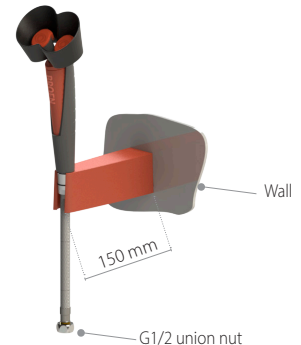

### Eye Shower with a Single 45° Head



No. 27 10 003 10 05 011

Weight: Approx. 0.5 kg

Table mounted eye shower with a single 45° head, incl. 1.5m inlet hose. Trigger operated including a built-in check valve and FLOWFIX for regulation of water flow (7 l/min). Min. working pressure: 1,5 bar.




### Eye Shower with a Single 45° Head



No. 27 10 003 10 03 011

Weight: Approx. 0.5 kg

Wall mounted eye shower with a single 45° head, incl. 1.5m inlet hose. Trigger operated including a built-in check valve and FLOWFIX for regulation of water flow (7 l/min). Min. working pressure: 1,5 bar.

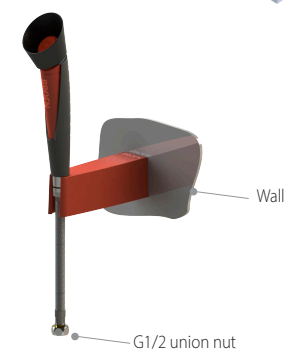

Other configurations are possible - please refer to ordering information on page 5.



## HAND-HELD EYE SHOWERS AND ACCESSORIES

### Eye Shower with a Single 45° Head

**No. 27 10 003 10 00 001**

Weight: Approx. 0.2 kg

Eye shower with a single head.  
Trigger operated including a built-in check valve and FLOWFIX for regulation of water flow (7 l/min).

Also available as straight model.



### Eye Shower with Two 45° Heads

**No. 27 10 009 10 00 001**

Weight: Approx. 0.2 kg

Eye shower with two angled heads.  
Trigger operated including a built-in check valve and FLOWFIX for regulation of water flow (14 l/min).

Also available as straight model.



## ACCESSORIES FOR HAND-HELD EYE SHOWERS

### Mounting Accessories (YY):

**No. 17 760 300**

Table base (M30x1.5, 60mm)



**No. 17 760 310**

30° Wall base



**No. 17 760 320**

Wall bracket



**No. 17 760 330**


Bracket for modular body shower



### Signs and Instructions (YY):

**No. 17 001 970**

Complete set with signs and instructions for eye shower

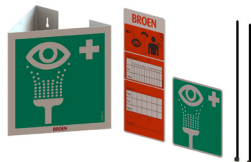

**No. 96M0339**

Sign for eye shower (sticker)



### Hoses and hose accessories (ZZ):

**No. 93H40561500**

1.5m Silicone hose with a Stainless Steel braiding and Swivel Joint for hand-held eye showers.



(other lengths on request).

**No. 17701760**

Counter weight/stop for hose



## ORDERING INFORMATION

For ordering a complete hand-held eye shower including accessories, please follow the following structure:

**27 10 XXX 10 YY ZZ 1**

#### Serial number

- 001 - Rinsing model, 45°
- 003 - Eye Shower, single, 45°
- 004 - Rinsing model, straight
- 006 - Eye shower, single, straight
- 009 - Eye shower, double, 45°
- 012 - Eye shower, double, straight

#### Country code

- 10 - International

#### \* Accessories

- 00 - None
- 01 - 17 001 970
- 02 - 96M0339
- 03 - 17 001 970 + Wall bracket
- 04 - 17 001 970 + 30° Wall Base
- 05 - 17 001 970 + Table base
- 06 - 96M0339 + Wall bracket
- 07 - 96M0339 + 30° Wall Base
- 08 - 96M0339 + Table base
- 09 - 96M0339 + Bracket for modular body shower

#### \* Hose length + Hose accessory

- 00 - None
- 01 - 1.5m with G1/2 union nut
- 02 - 1.5m with G1/2 union nut + stop
- 03 - 2m with G1/2 union nut
- 04 - 2m with G1/2 union nut + stop

\* Please request document No. 99G0045 for a complete list of available hoses and accessories.

For more information and technical datasheets please contact your BROEN-LAB sales representative.

## EYE/FACE SHOWERS

BROEN-LAB has developed the Eye/Face showers, which comply to EN15154-2 and ANSI Z358.1 norms and also offer many improved features.



### MATERIALS WITH MEDIA CONTACT

Plastics: TPE  
Metals: Brass  
Rubber: EPDM

## SAFETY

- Flow rate in according to EN15154-2
- BROEN-LAB recommends the installation of a Thermostatic Mixing Valve (TMV) with safety shower equipment for safe operation, see page 11
- Sign according to ISO 3864-1


## OPERATION

- Stay-open dust cap
- Plug and play assembly
- Spray pattern according to EN15154-2
- Stronger wall flange with 6 holes for wall mounted models

## RELIABILITY

- Easy-to-clean spray head
- Anti-limescale design and material
- Improved functionality
- Long lifetime

### Eye/Face Shower with Bowl, Wall Mounted

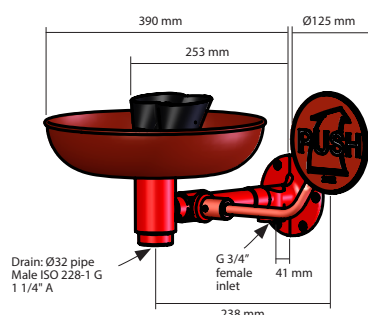

Weight: Approx. 4.5 kg

Eye/Face shower with two heads in brass with chemical resistant BROEN-LAB Polycoat and Stainless Steel bowl. Wall mounted for exposed pipe work. The water amount is regulated by a built-in FLOWFIX (14 l/min). Min. working pressure: 1 bar.

#### No. 17 528 009

No. 17 528 009E  
Includes a 700 mm connection hose with union nut 1/2\" inlet. More information, see page 13

#### Other Models

No. 17 533 009

Male inlet



#### Anti-corrosive

No. 17 528 004

G 1/2\" female inlet




No. 17 533 004

Male inlet



### Eye Shower, Wall Mounted

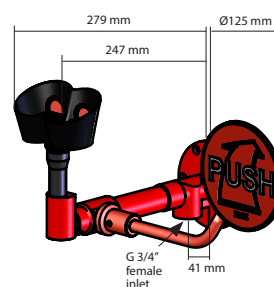


Weight: Approx. 4.5 kg

Eye shower, wall mounted in brass with chemical resistant BROEN-LAB Polycoat. The water is regulated with a built-in FLOWFIX (14 l/min). Min. working pressure: 1 bar.

#### No. 17 550 009

No. 17 550 009E  
Includes a 700 mm connection hose with union nut 1/2\" inlet. More information, see page 13

### Eye Shower, Table Mounted, Swing Down

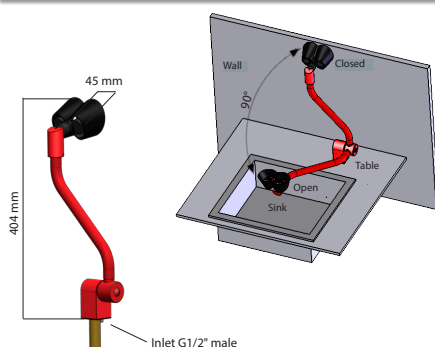

Weight: Approx. 4.5 kg


Eye shower, table mounted, swing down in brass with chemical resistant BROEN-LAB Polycoat for mounting on back side of the sink. The water amount is regulated with built-in FLOWFIX (26 l/min). Min. working pressure: 2 bar.

#### No. 17 527 009


#### Other Models

No. 17 542 009

For swivel wall mounting

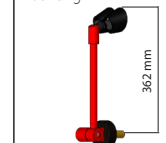

### Eye Shower, Table Mounted, Swivel Joint

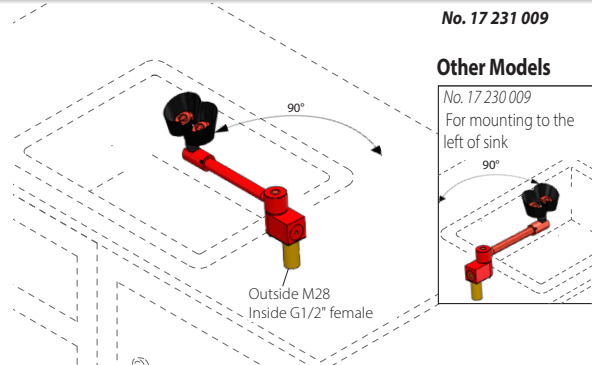

Weight: Approx. 4 kg


Eye shower, table mounted, swivel joint in brass with chemical resistant BROEN-LAB Polycoat for mounting on the right side of sink. The water amount is regulated with the built-in FLOWFIX (26 l/min). Min. working pressure: 1.5 bar.

#### No. 17 231 009

#### Other Models

No. 17 230 009

For mounting to the left of sink

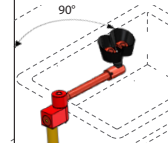

## SPECIAL FEATURES

### GENERAL:


- Choose between BROEN-LAB Polycoated surface or Anti Corrosive.
- Cabinet and cabinet door in Stainless Steel AISI 304.
- Disability-Friendly models comply with ANSI A117.1-2003.
- Easy to mount and install.
- Cabinet can also be surface mounted (except 17 155 009).

Choose among different body shower outlets by ordering recessed units and outlets separately to fit your exact needs of installation and application.

E.g. : 17 154 009 + 17 110 009

### BODY SHOWER:

- Strong integrated handle for operating body shower.
- Free choice of combination between ceiling or wall mounted body shower outlet.
- Premounted valve for body shower.

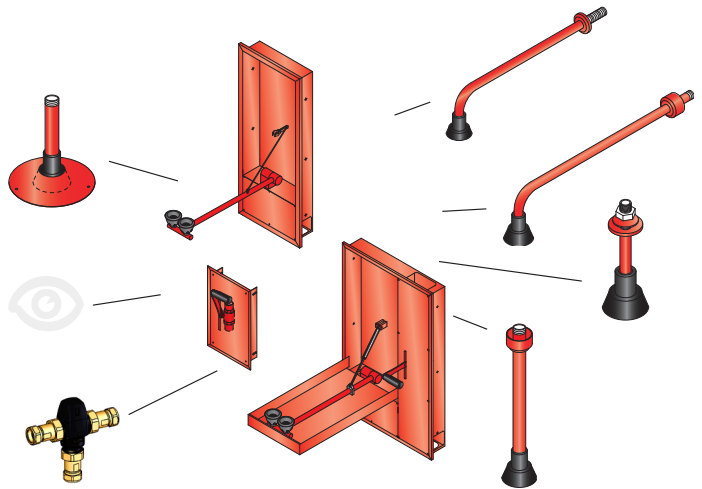


### EYE WASH:

- Cabinet and cabinet door in Stainless Steel AISI 304.
- Integrated drain (on models with door).
- Cabinet door catches the water and prevents it from splashing on the floor.
- Filter BALLOFIX® ensures clean shower heads.

### Body Shower Outlet, Ceiling Mounted, Concealed




No. 17 100 009

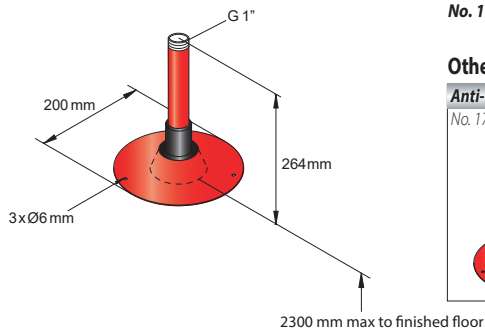

#### Other Models

Anti-corrosive

No. 17 100 004

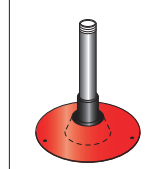

Weight: Approx. 1 kg


Body shower outlet, ceiling mounted, concealed in Stainless Steel with self-draining shower head. Inlet G 1" male.

Min. working pressure: 1.5 bar.

### Body Shower Outlet, Ceiling Mounted



No. 17 110 009

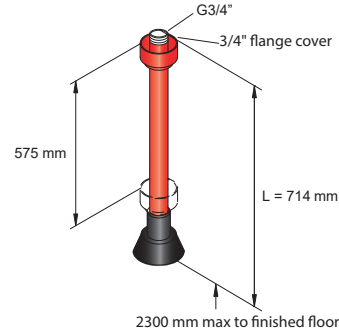

#### Other Models

No. 17 089 009

Male inlet and fixed flange

G 3/4" male

L = 271 mm



#### Anti-corrosive

Same as ill. Polycoat but anti-corrosive. No. 17 110 004

Male inlet and fixed flange No. 17 089 004

Weight: Approx. 1.5 kg


Body shower outlet, ceiling mounted, in brass with chemical resistant BROEN-LAB Polycoat with self-draining shower head. Fits different ceiling heights by adjusting the flange. Inlet G 3/4" male.

Min. working pressure: 1.5 bar.

### Body Shower Outlet, Wall Mounted



No. 17 120 009

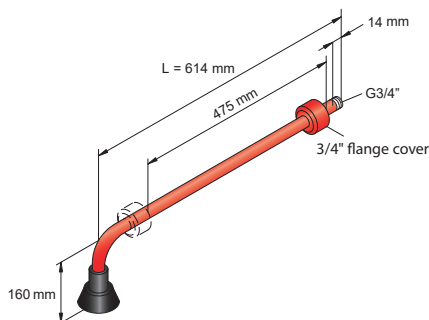


#### Other Models

No. 17 088 009

Male inlet and fixed flange

G 3/4" male

L = 510 mm

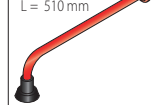

Weight: Approx. 1.5 kg

Body shower outlet, wall mounted, in brass with chemical resistant BROEN-LAB Polycoat with self-draining shower head. Inlet G 3/4" male. With adjustable flange.

Min. working pressure: 1.5 bar.


### Body Shower Outlet, Wall Mounted

Anti-corrosive



### Disability-Friendly

No. 17 130 004

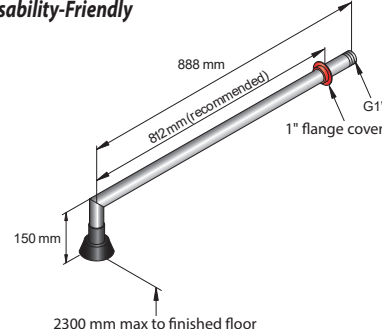

#### Other Models

No. 17 088 004

Male inlet and fixed flange. Stanniloy

L = 510 mm

G 3/4" male



Weight: Approx. 1.5 kg

Wall mounted body shower outlet in Stainless Steel AISI 304, with self-draining shower head. Inlet G 1" male. With adjustable flange.

Dimensions comply with the requirements of ANSI A117.1-2003.

Min. working pressure: 1.5 bar.



# RECESSED EMERGENCY SHOWERS

## Eye Shower, Wall Recessed




No. 17 150 009

### Other Models

#### Anti-corrosive

Same as ill. Polycoat but anti-corrosive cabinet.  
No. 17 160 009

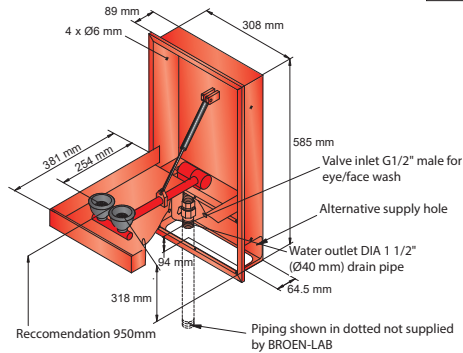

Weight: Approx. 8 kg

Eye shower, wall recessed with cabinet with chemical resistant BROEN-LAB Polycoat, two outlet heads and built-in FLOWFIX (14 l/min). Cabinet door activates the eye wash and catches the water. Min. working pressure: 2 bar.

## Eye Shower, Wall Recessed




No. 17 152 009

### Disability-Friendly

### Other Models

#### Anti-corrosive

Same as ill. Polycoat but anti-corrosive cabinet.  
No. 17 162 009

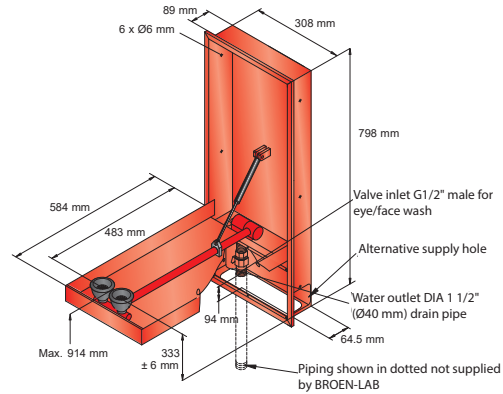

Weight: Approx. 10 kg

Eye shower, wall recessed, disability-friendly with cabinet with chemical resistant BROEN-LAB Polycoat, two outlet heads and built-in FLOWFIX (14 l/min). Cabinet door activates the eye wash and catches the water. Min. working pressure: 2 bar.

## Eye and Body Shower, Wall Recessed




No. 17 154 009

### Disability-Friendly

### Other Models

#### Anti-corrosive

Same as ill. Polycoat but anti-corrosive cabinet.  
No. 17 164 009

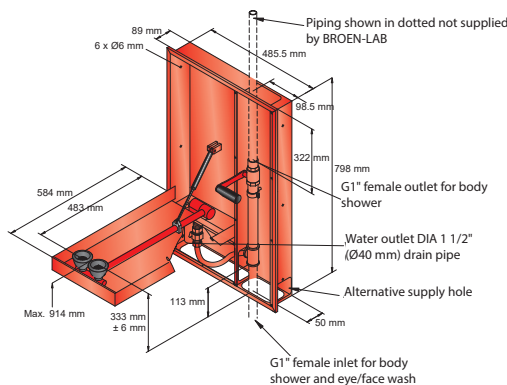

Weight: Approx. 16 kg

Eye shower, with 1" valve for body shower, wall recessed, disability-friendly with cabinet with chemical resistant BROEN-LAB Polycoat, two outlet heads and built-in FLOWFIX (14 l/min). Cabinet door activates the eye wash and catches the water. Integrated hand activating 1" valve for body shower with various options for shower outlet. Min. working pressure: 2 bar.

## Body Shower, Wall Recessed, Hand Activator




No. 17 155 009

### Other Models

#### Anti-corrosive

Same as ill. Polycoat but anti-corrosive cabinet.  
No. 17 165 009

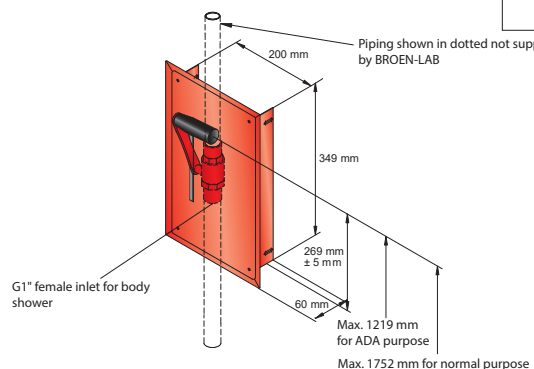


Weight: Approx. 2.5 kg

Wall recessed cabinet with chemical resistant BROEN-LAB Polycoat and hand activator for body shower with various options for shower outlet. 1" valve for 3/4" or 1" pipe. Min. working pressure: 1 bar.

## BODY SHOWERS

BROEN-LAB has developed the body showers, which comply to EN15154-1 and ANSI Z358.1 norms and also offer many improved features.

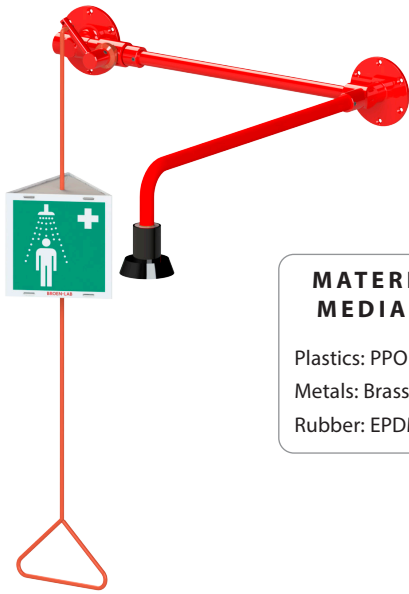

### MATERIALS WITH MEDIA CONTACT

Plastics: PPO  
Metals: Brass  
Rubber: EPDM

## SAFETY

- Flow rate in according to EN15154-1
- Spray pattern according to the EN 15154-1
- Sign according to ISO 3864-1
- BROEB-LAB recommends the installation of a Thermo-static Mixing Valve (TMV) with safety shower equipment for safe operation, see page 11

## OPERATION


- Plug and Play assembly
- Stronger wall flange with 6 holes for wall mounted models with exposed piping

## RELIABILITY

- Plug and Play Spare Parts
- Anti-limescale design and material
- Improved functionality
- Long lifetime

### Body Shower, Ceiling Mounted


No. 17 580 009

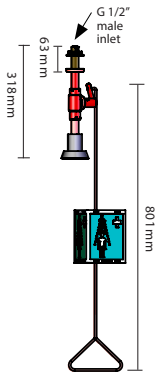

Weight: Approx. 2 kg  
Body shower, ceiling mounted with self-draining head. In brass with chemical resistant BROEN-LAB Polycoat.  
Min. working pressure: 1 bar.

### Body Shower, Wall Mounted

No. 17 581 009

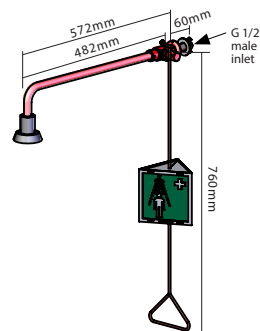


Supporting wall bracket can be used for mounting (see page 11)

Weight: Approx. 3 kg  
Body shower, wall mounted with self-draining shower head. In brass with chemical resistant BROEN-LAB Polycoat.  
Min. working pressure: 1 bar.

### Body Shower for Exposed Pipework, Wall Mounted

No. 17 582 009

No. 17 582 009E  
Includes a 700 mm connection hose with union nut 1/2" inlet. More information, see page 13

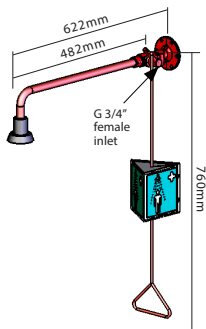


Weight: Approx. 3 kg  
Body shower for exposed pipe-work with self-draining head. In brass with chemical resistant BROEN-LAB Polycoat.  
Min. working pressure: 1 bar.

Supporting wall flange can be used for mounting (see page 11)

### Body Shower, Wall Mounted above Door

No. 17 587 009

No. 17 587 009E  
Includes a 700 mm connection hose with union nut 1/2" inlet. More information, see page 13

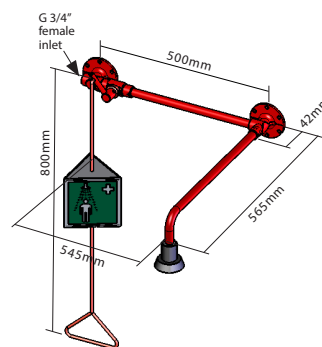


Weight: Approx. 4.5 kg  
Body shower, wall mounted above door for exposed pipework with self-draining head. In brass with chemical resistant BROEN-LAB Polycoat. Valve to be mounted left or right of door.  
Min. working pressure: 1 bar.

## COMBINATION SHOWERS

### Body Shower with Eye/Face Shower and Bowl, Modular

#### No. 17656009

No. 17656009E  
Includes a 700 mm connection hose with union nut 1/2" inlet.  
More information, see page 13

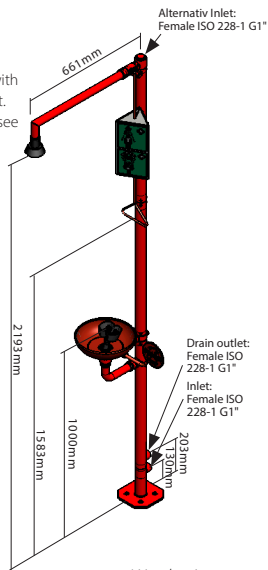


Weight: Approx. 17 kg

Body shower with eye/face shower and bowl, modular. In brass and SS with chemical resistant red BROEN-LAB Polycoat with self-draining shower head. Eye wash with built-in FLOWFIX for regulation of water flow (14 l/min). Bowl in Stainless Steel with integrated outlet.  
Min. working pressure: 1.5 bar.

Anti-corrosive - Modular

No. 17656004

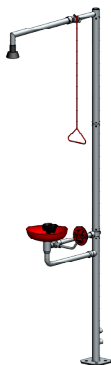


Same as ill. Polycoat but anti-corrosive.

### Body Shower with Eye/Face Shower and Bowl, Modular

#### Disability-Friendly

No. 17657009

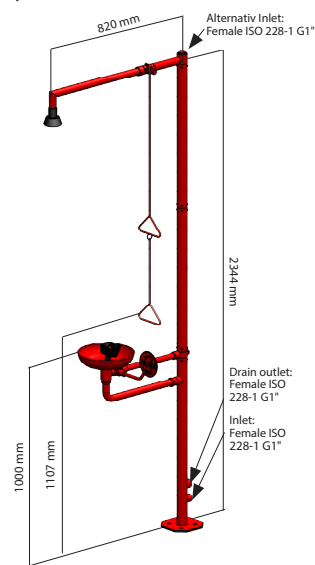

Weight: Approx. 22 kg

Body shower with eye/face shower and bowl, modular. In brass and SS with chemical resistant red BROEN-LAB Polycoat with self-draining shower head. Eye wash with built-in FLOWFIX for regulation of water flow (26 l/min). Bowl in Stainless Steel with integrated outlet.  
Min. working pressure: 1.5 bar.

Anti-corrosive - Modular

#### Disability-Friendly

No. 17657004




Same as ill. Polycoat but anti-corrosive.

### Body Shower with Eye/Face Shower and Bowl, Modular

#### No. 17556009


No. 17556009E  
Includes a 700 mm connection hose with union nut 1/2" inlet.  
More information, see page 13

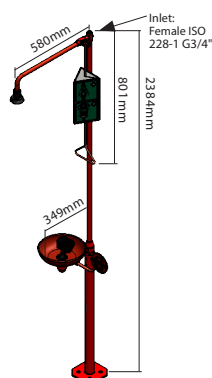

Weight: Approx. 17 kg

Body shower with eye/face shower and bowl, modular. In brass and SS with chemical resistant red BROEN-LAB Polycoat with self-draining shower head. Eye wash with built-in FLOWFIX for regulation of water flow (14 l/min). Bowl in Stainless Steel.  
Min. working pressure: 1.5 bar.


No. 17956609  
With hand-held eye shower

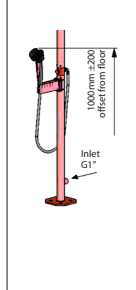

### Body Shower with Eye/Face Shower and Bowl, Modular, Wall Mounted

#### No. 17552009



No. 17552009E  
Includes a 700 mm connection hose with union nut 1/2" inlet.  
More information, see page 13

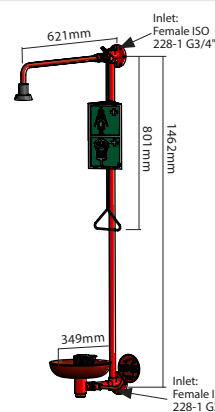

Weight: Approx. 8 kg

Body shower with eye/face shower and bowl, modular. In brass and SS with chemical resistant red BROEN-LAB Polycoat with self-draining shower head. Eye wash with built-in FLOWFIX for regulation of water flow (14 l/min). Bowl in Stainless Steel.  
Min. working pressure: 1.5 bar.

#### Other Models

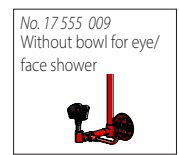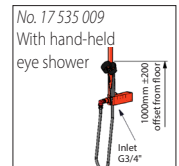



## Thermostatic Mixing Valves for Emergency Shower Systems

**BROEN-LAB TMVs for eye and body showers are thermostatic mixing valves designed for delivering tepid water (20-25°C) for eye/hand or body showers.**


The solution is designed for the operation of one eye/hand shower respectively one body shower.

The mixed temperature is pre set to tepid water from factory.

**The TMV incorporates in one body the following features:**

- Fail-safe on hot side: If the cold water is interrupted the built-in anti-scalding device will automatically cease the flow of hot water. **AUTOMATIC OPERATION**
- Fail-safe on cold side: If the hot water is interrupted the built-in intentional leak will allow cold water to flow at a sufficient flowrate to ensure operation. **THIS WORKS AUTOMATICALLY AS WELL**
- Including check valves on both hot and cold supply

## Thermostatic Mixing Valves for Eye Shower

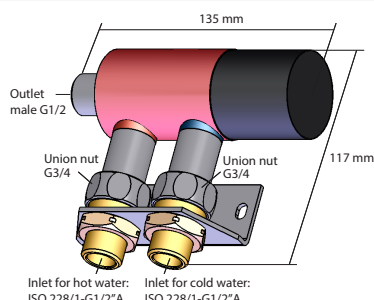


Weight: Approx. 1.5 kg

**No. 17140039:** Model with G1/2" male inlet  
Thermostatic mixing valve for eye shower with bracket.  
Min. working pressure: 1 bar.

**No. 17 140 009:** Model with G3/4" male inlet  
Thermostatic mixing valve for eye shower without bracket.  
Min. working pressure: 1 bar.

## Thermostatic Mixing Valves for Body- and Combination Showers

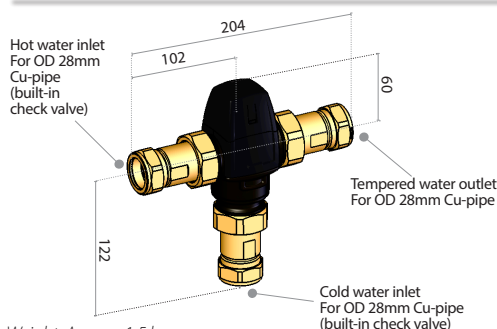

Weight: Approx. 1.5 kg

**No. 17 512 009:**  
Thermostatic mixing valve for body and Eye shower.  
Min. working pressure: 1.5 bar.

**No. 17 513 009:** Model with R1" inlet/outlet connections

## Signs for Emergency Shower Systems

### Nr. 17 001 920

(Weight: Approx. 0.23 kg)

Sign (pull rod mounted) and instruction card for Body and Eye Shower



### Nr. 17 002 27

(Weight: Approx. 0.03 kg)

Instruction and Service Sign - Body Shower



### Nr. 96M0042

(Weight: Approx. 0.01 kg)

Sign for body shower (sticker)



### Nr. 17 001 910

(Weight: Approx. 0.13 kg)

Sign (pull rod mounted) and instruction card for Body Shower



### Nr. 17 002 28

(Weight: Approx. 0.03 kg)

Instruction and Service Sign - Body and Eye Shower



### Nr. 96M0043

(Weight: Approx. 0.01 kg)

Sign for eye shower (sticker)



### Nr. 17 002 010

(Weight: Approx. 0.15 kg)

Sign (wall mounted) and instruction card for Body and Eye Shower



### Nr. 17 002 29

(Weight: Approx. 0.03 kg)

Instruction and Service Sign - Eye Shower fountain



### Nr. 96M0345

(Weight: Approx. 0.01 kg)

Sign for direction (sticker)



### Nr. 17 001 900

(Weight: Approx. 0.08 kg)

Sign (wall mounted) and instruction card for Eye Shower



### Nr. 17 002 30

(Weight: Approx. 0.03 kg)

Instruction and Service Sign - Eye Shower - Hand-Held



### Nr. 17 002 000

(Weight: Approx. 0.08 kg)

Sign (wall mounted) and instruction card for Body Shower



### Nr. 17 002 31


(Weight: Approx. 0.03 kg)

Instruction and Service Sign - Body Shower with Hand-Held Eye Shower






### Cable

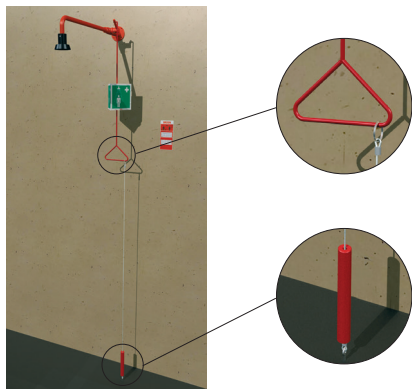

#### No. 17 901 64

A cable may be installed as a supplement to the pull rod solution i.e. the cable enables the user to activate the shower at floor level.


### Wall bracket for body showers

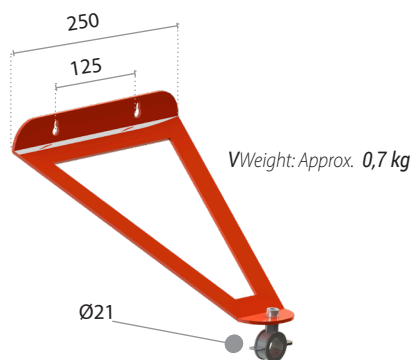

#### Nr. 17 001 079

Wall bracket for support of wall mounted body showers (size 1/2"). Other sizes on request

### Foot Panel for Operation of Free Standing Body Shower

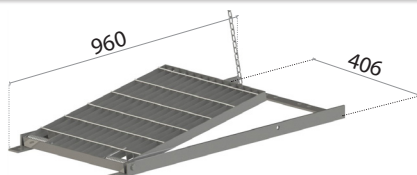


Weight: Approx. 23 kg

#### No. 17 919 009 - Stainless Steel

#### No. 17 834 009 - Galvanized Steel

Foot panel for operation of body showers. Possible to combine with 17 909.909 if mounted 180° opposite. Delivered with all mounting accessories

### Foot Treadle for Operation of Eye Showers

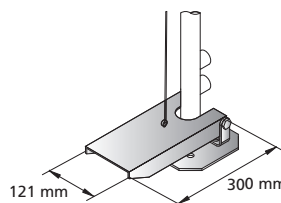

Weight: Approx. 1 kg

#### No. 17 909 909

Foot treadle in aluminium and SS for operation of eye showers. Possible to combine with 17 919.009 if mounted 180° opposite. Delivered with all mounting accessories.

### Alarm for Emergency Shower



#### No. 17 851 95

Signal and light alarm for emergency showers.

- 112dB(A) output
- IP66 protection
- 32 alarm tones

Working temperature: -20°C to 70°C

Power supply: 230VAC, 50/60Hz

For operation in Non Ex Proof areas (also available for Ex Proof areas)

### Flow switch



Weight: Approx. 0,5 kg

- Good reliability
- Low pressure loss
- Dirt resistant

Item no.	Material	Size*	PN, bar
96P0062	brass	G3/4	25
96P0063	brass	G 1 1/4	25
96P0066	SS	G3/4	25
96P0068	SS	G 1 1/4	25

\* Please note that size of the flow switch should be a size bigger than the inlet of the emergency shower

- Hermetic separation of electrical and hydraulic components

When there is a stream of water in the connection piping to the emergency shower the paddle in the flow switch is triggering an electrical reed switch. This switch can be connected to any type of warning system, such as an alarm central, a horn, a light or similar.

For operation in Non Ex Proof areas (also available for Ex Proof areas)

For more information and technical datasheets please contact your BROEN-LAB sales representative.



The standard 700 mm connection hose with G1/2" union nut is suitable for all new models with female G3/4". For other requirements, do not hesitate to contact BROEN-LAB.



## Easy to install on current pipes

At BROEN-LAB we focus on quality in every detail to help you work smarter and reduce installation time.



## No teflon tape or hemp is needed

To reduce the installation time the connection is o-ring sealed. Therefore no Teflon tape is needed.




## O-ring sealing

Additionally the O-ring connection secure less risk of leakage and raise the level of safety.



See how easy to connect


## Connection hoses for Emergency Showers

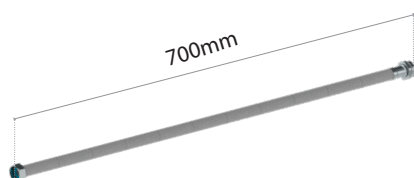

- No. 93i59370700 (Male G1" - Female G1" union nut connection)
- No. 93i42370700 (Male G1" - Female G3/4" union nut connection)
- No. 93i42330700 (Male G3/4" - Female G3/4" union nut connection)
- No. 93i42370700 (Male G3/4" - Female G1" union nut connection)

Material: Silicone  
Braiding: Stainless Steel  
Inside diameter: Ø15mm  
Outside diameter: Ø21,5mm  
Max.working pressure: 10 bar  
Max.working temperature: 90°C

## Drain kit for Emergency Showers

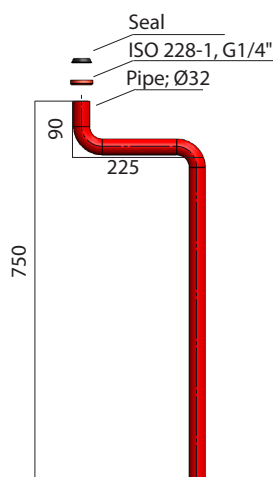

No. 17 556 579

For extra drainage under bowl on body showers with eye/face wash and bowl.

Compatible with:  
17 552 009  
17 528 099  
17 541 009  
17 556 009  
17 533 009



# TEST EQUIPMENT FOR EMERGENCY SHOWER SYSTEMS

## Test Trolley for Emergency Showers

### No. 17 851 929, model with integrated pump

\*With 8 mm connection cable and international plug


Test trolley with integrated pump.

Drainage hose is included. Capacity: 90 l

Can be used for heated and unheated emergency showers.

Weight: Approx. 22 kg

For operation in Non Ex Proof areas.

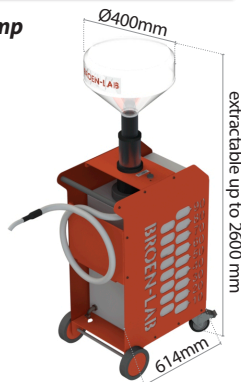

## Test Funnel for Emergency Showers

### No. 17 001 08

Test funnel for testing of emergency showers in PC (polycarbonate) with 2.2m hose.

### No. 27 99 008 10 01

Measuring Cylinder for checking Eye Wash Stations



## FIRST AID KIT

**P7330051000002**



Standard First Aid Kit which contains

- Clean pads
- Cool pads
- Stop dressing
- Elastic band aids x 2
- Sterile eye wash
- pH eye wash
- Signs
- instruction sheets
- First Aid Plans

It's also possible to configure you own First Aid Kit using the refill components beneath. For help please contact your BROEN-LAB representative.

## REFILLS

### Quick Clean

**P7530071000000**



Sterile wipes are ideal to clean skin and wound prior to dressing.

### Quick Stop

**P7530081000000**



Wound dressing kit used to protect wound form further harm.

### Quick Cool

**P7530061000000**



Burn gel that eases the pain if the skin has been in contact with hot surfaces or substances.

### pH Neutral Eye Wash

**P7170031000000**



PH Neutral features phosphate buffer property with documented efficiency in chemicals-based eye injuries.

### Elastic Plaster

**P7530061000000**



Elastic plasters that conform the movement of the skin so they are comfortable to use. The plasters are extremely skin-friendly and allow the skin to breathe.

### Sterile Eye Wash

**P7170021000000**



Sterile Eye wash for optimum assistance if dirt, dust or other substances come in contact with eyes.

For more information and technical datasheets please contact your BROEN-LAB sales representative.

## INTRODUCTION

**Thousands of industrial accidents happen every year in spite of precautionary safety measures, accidents that are caused by human exposure to heat (fire, boiling water and oil, steam, etc.), to corrosive agents (fluid acids and caustics) and/or to influence by dust and corrosive aerosols (metal, dirt, fluids).**



These accidents cause the injured person acute pain, disablement and may even lead to loss of life. However, experience shows that the consequences and human suffering following such accidents can be considerably limited if water is applied directly to the affected areas instantly after the exposure.

It is a fact of life that the risk lurks in many places: in the laboratory, in the kitchen, in the warehouse, in the chemical plant, in the workshop or generally speaking in any place where heat, corrosion or dust can cause serious human injury. For these reasons we strongly recommend the installation of adequate emergency shower systems as a means of limiting these risks.

## Location and Operation of Emergency Showers

As a basic rule the shower must be located as close to the workplace as possible so that the general first aid instructions of flushing the affected area immediately after the accident has occurred can be followed.

Please note that if strong acids or caustics are involved immediate action is demanded. When installing the emergency shower system please also be aware that it must be easy to locate and access at all times for potential users.

BROEN-LAB Emergency Showers can be recognized by their red and white stripes on body and eye showers, red pull rods and luminescent signs and symbols for the clear identification of emergency shower.


*Mobile and portable emergency showers - see page 18*



*Heated body showers - see page 16-17*

## Selecting the Right Shower

There are many different conditions that determine the type of emergency shower to be installed. It is important that the chosen unit is designed and tested according to the existing standard as well as suits the specific application and location - **and environment!**

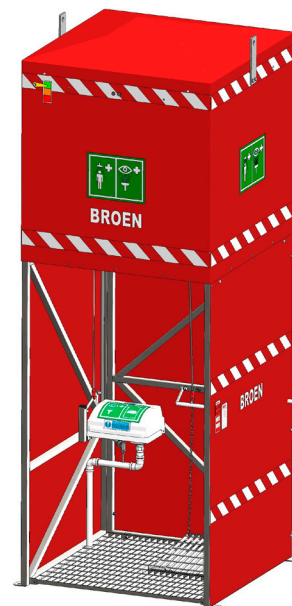

*Tank showers - see page 19*

**In this section you can find a wide range of BROEN-LAB Emergency Showers that can be installed outside and be used in areas where outside temperatures go as low as -35°C.**




## HEATED COMBINATION SHOWERS

Frost protected outdoor combination safety showers, trace tape heated and pre-insulated. Delivered with eye/facewash fountain and additional eyebath diffuser on flexible hose to treat unconscious or disabled personnel.

Available in Galvanized or Stainless Steel.


### Heated Body Shower - Free Standing with Eye Wash Fountain

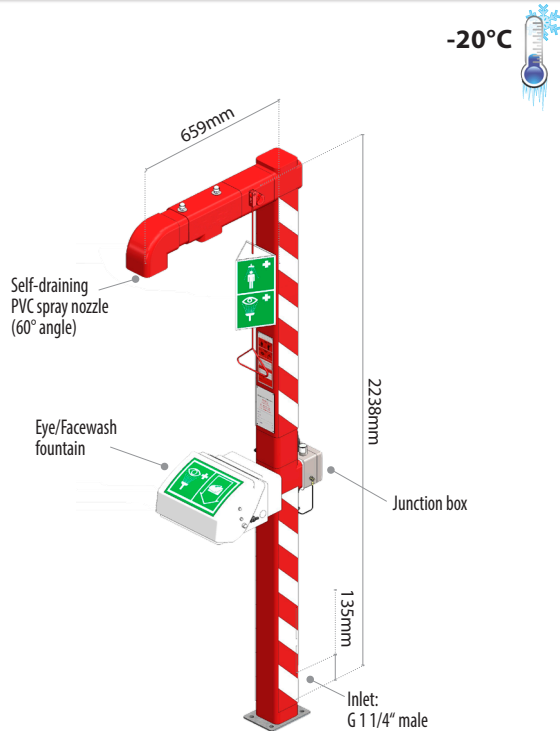


#### TECHNICAL SPECIFICATIONS

Min. working pressure: 2 Bar  
Min. working temperature: -20°C  
Water Flow: 115 l/min.  
Supply voltage: 230 VAC, 50/60 Hz

*Heated underground extension pipe*  
(ordered separately, see page 20 for more information)

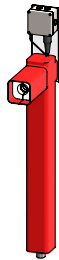

### Heated Body Shower - Free Standing with Eye Wash Fountain

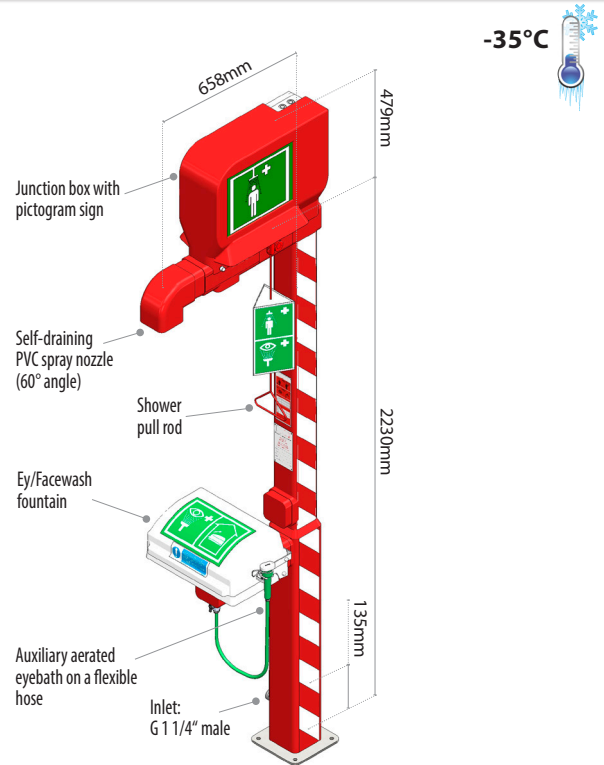


#### TECHNICAL SPECIFICATIONS

Min. working pressure: 2 Bar  
Min. working temperature: -35°C  
Water Flow: 115 l/min.  
Supply voltage: 230VAC, 50/60Hz

*Heated underground extension pipe*  
(ordered separately, see page 20 for more information)



## ORDERING INFORMATION

BROEN-LAB item no.	Material	 Ex-proof*	Min. working temperature	Eye wash fountain	Mounting	Heated underground extension pipe
H7 41 002 10 00 321	Galvanized Steel	Yes	-20°C	Yes	Floor	Available
H7 41 002 10 00 221	Stainless Steel	Yes	-20°C	Yes	Floor	Available
H7 41 001 10 00 341	Galvanized Steel	Yes	-35°C	Yes	Floor	Available
H7 41 001 10 00 241	Stainless Steel	Yes	-35°C	Yes	Floor	Available
H7 41 001 10 00 331	Galvanized Steel	No	-35°C	Yes	Floor	Available
H7 41 001 10 00 231	Stainless Steel	No	-35°C	Yes	Floor	Available

\*Ex-proof zone 1 and zone 2.


For more information and technical datasheets please contact your BROEN-LAB sales representative.



## HEATED BODY SHOWERS

Frost protected outdoor combination safety showers, trace tape heated and pre-insulated.  
Available in Galvanized or Stainless Steel.


### Heated Body Shower - Free Standing

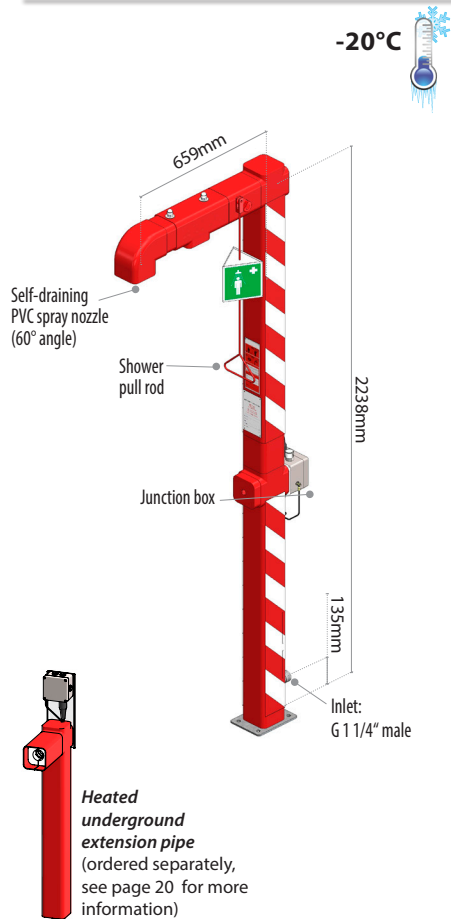

#### TECHNICAL SPECIFICATIONS

Min. working pressure: 2 Bar  
Min. working temperature: -20°C  
Water Flow: 115 l/min.  
Supply voltage: 230VAC, 50/60Hz


### Heated Body Shower - Free Standing

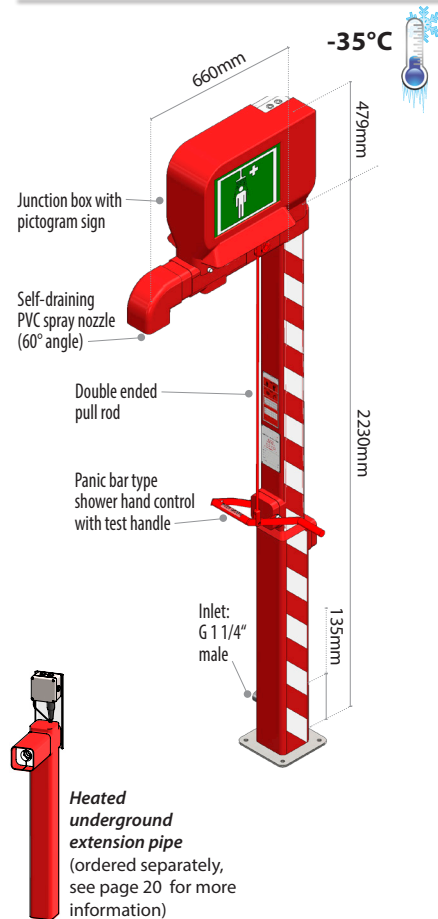

#### TECHNICAL SPECIFICATIONS

Min. working pressure: 2 Bar  
Min. working temperature: -35°C  
Water Flow: 115 l/min.  
Supply voltage: 230VAC, 50/60Hz

### Heated Body Shower - Wall Mounted

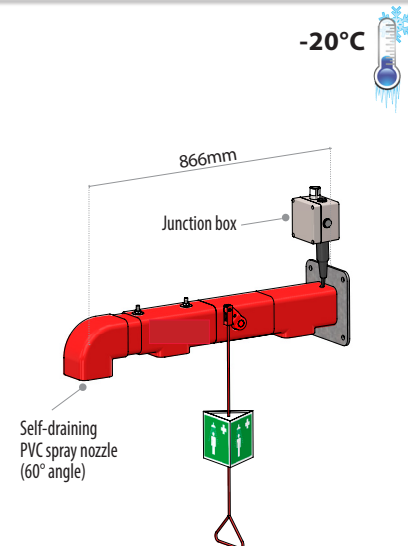

#### TECHNICAL SPECIFICATIONS

Min. working pressure: 2 Bar  
Min. working temperature: -20°C  
Water Flow: 115 l/min.  
Supply voltage: 230VAC, 50/60Hz

## ORDERING INFORMATION


BROEN-LAB item no.	Material	Ex-proof*	Min. working temperature	Eye wash fountain	Mounting	Heated underground extension pipe
H7 21 002 10 00 321	Galvanized Steel	Yes	-20°C	No	Floor	Available
H7 21 002 10 00 221	Stainless Steel	Yes	-20°C	No	Floor	Available
H7 21 001 10 00 341	Galvanized Steel	Yes	-35°C	No	Floor	Available
H7 21 001 10 00 241	Stainless Steel	Yes	-35°C	No	Floor	Available
H7 21 001 10 00 331	Galvanized Steel	No	-35°C	No	Floor	Available
H7 21 001 10 00 231	Stainless Steel	No	-35°C	No	Floor	Available
H7 23 001 10 00 321	Galvanized Steel	Yes	-20°C	-	Wall	-
H7 23 001 10 00 221	Stainless Steel	Yes	-20°C	-	Wall	-

\*Ex-proof zone 1 and zone 2.

For more information and technical datasheets please contact your BROEN-LAB sales representative.

## HEATED EYE/FACE SHOWERS

### Heated Eye Wash Fountain, Pedestal Mounted

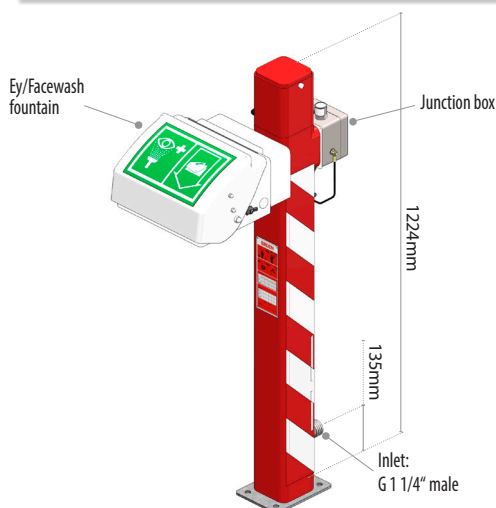

Outdoor pedestal mounted emergency eye/facewash fountain, trace tape heated and pre-insulated. Delivered with “unbreakable” integral lid, aerated diffusers, flow regulator and self-adhesive sign.

Available in Galvanized or Stainless Steel.

#### TECHNICAL SPECIFICATIONS

Min. working pressure: 2 Bar  
Min. working temperature: -20°C  
Water Flow: 20 l/min.  
Supply voltage: 230VAC, 50/60Hz

Heated underground extension pipe (ordered separately, see page 20 for more information)



### ORDERING INFORMATION


BROEN-LAB item no.	Material	Ex-proof*	Min. working temperature	Mounting	Heated underground extension pipe
H7 31 001 10 00 321	Galvanized Steel	Yes	-20°C	Floor	Available
H7 31 001 10 00 221	Stainless Steel	Yes	-20°C	Floor	Available

\*Ex-proof zone 1 and zone 2.

## MOBILE & PORTABLE EMERGENCY SHOWERS

Mobile Emergency Safety Shower is designed as a quick response, first-aid device for use where a plumbed-in unit is not available and must therefore be backed up by a method of further decontamination.

### Mobile Self Contained Emergency Safety Shower

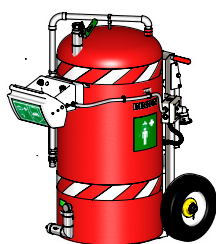

#### TECHNICAL SPECIFICATIONS

Capacity: 114 l  
Min. working temperature: -10°C  
Water Flow: 40 l/min.  
Supply voltage: 230VAC, 50/60Hz

#### ORDERING INFORMATION

BROEN-LAB item no.	Material	Ex-proof*	Min. working temperature
H7 47 002 10 00 311	Galvanized & SS	No	-10°C
H7 47 002 10 00 321	Galvanized & SS	Yes	-10°C
H7 47 002 10 00 301	Galvanized & SS	No	+1°C (unheated)

\*Ex-proof zone 1 and zone 2.

### Portable Emergency Eye, Face and Body Shower (Unheated)




No. H7170011000303

#### TECHNICAL SPECIFICATIONS

Capacity: 15 l  
Min. working temperature: +1°C  
Water Flow: 8 l/min.  
Material: Stainless Steel cylinder

### Mobile Self Contained 1200 L Shower and Eyewash incl. Pump

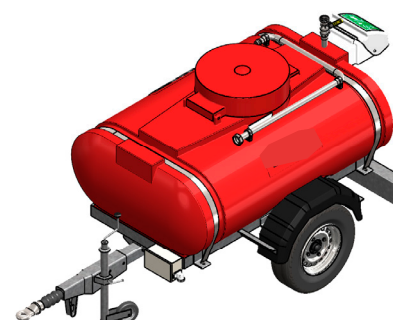

No. H7470031000331

#### TECHNICAL SPECIFICATIONS

Capacity: 1200 l  
Min. working temperature: -10°C  
Water Flow: 75 l/min.  
Supply voltage: 230VAC, 50/60Hz  
Material:  
- Galvanized mild steel trailer  
- Red polyethylene tank

## HEATED TANK SHOWERS


A range of heavy duty emergency tank showers, for use when warm water is required, or a constant supply of water cannot be guaranteed.

These units are essential where there is a problem with insufficient water pressure or inadequate flow rate. This improved range incorporates models to suit a wide range of situations.

### Body Shower with Eye Wash Fountain & 350 l tank



### Body Shower with Eye Wash Fountain & 1200 l tank

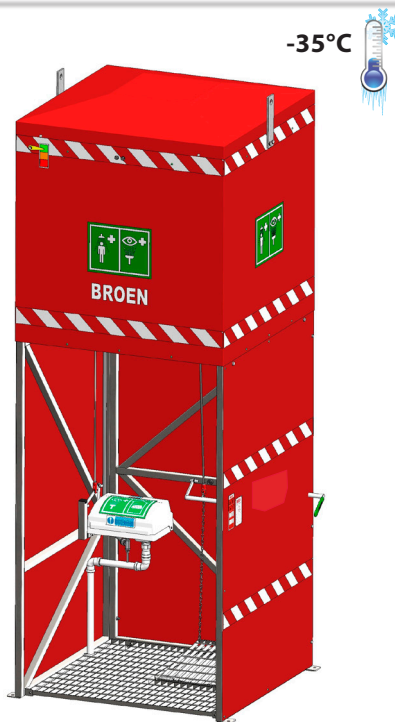

#### Standard Configuration

- 350 or 1200 l tank in medium density polythene
- Support legs and frame in galvanized or Stainless Steel
- Rear and side GRP panelling
- Visual level indicator
- Body shower
- Eye/facewash fountain
- Shower control panic bar
- Foot panel for body shower

#### TECHNICAL SPECIFICATIONS

Min. working temperature: -35°C  
Water Flow: 75 l/min.  
Supply voltage: 230VAC, 50/60Hz

*\*750l tank shower is available on request*

#### Tank showers are available in three sizes:

The **350 liter** model supplies a guaranteed flow of warm water for several minutes.

With the mains supply permanently connected, the **750\* liter** model will deliver warm water for 15 minutes, recommended by ANSI Z 358.1 and EN15154-1&2 standards.


The **1200 liter** models will meet the mentioned standards even without an external supply connected.

Available in Galvanized or Stainless Steel.

For areas with high outside temperatures models with chillers can be supplied.

Unheated models are also available.

## ORDERING INFORMATION



BROEN-LAB item no.	Material	 Ex-proof*	Min. working temperature	Eye wash fountain	Tank volume
H7 61 001 10 01 341	Galvanized Steel	Yes	-35°C	Yes	350 l
H7 61 001 10 01 241	Stainless Steel	Yes	-35°C	Yes	350 l
H7 61 001 10 01 331	Galvanized Steel	No	-35°C	Yes	350 l
H7 61 001 10 01 231	Stainless Steel	No	-35°C	Yes	350 l
H7 61 003 10 01 341	Galvanized Steel	Yes	-35°C	Yes	1200 l
H7 61 003 10 01 241	Stainless Steel	Yes	-35°C	Yes	1200 l
H7 61 003 10 01 331	Galvanized Steel	No	-35°C	Yes	1200 l
H7 61 003 10 01 231	Stainless Steel	No	-35°C	Yes	1200 l

\*Ex-proof zone 1 and zone 2.

For more information and technical datasheets please contact your BROEN-LAB sales representative.

## ACCESSORIES FOR HEATED SHOWERS

### Heated Underground Extension Pipe for Heated Body Showers

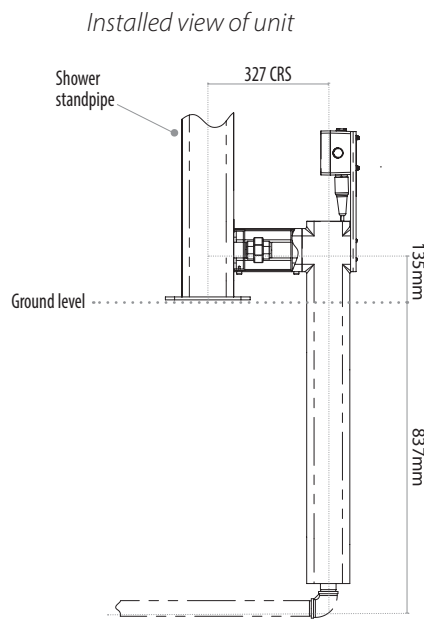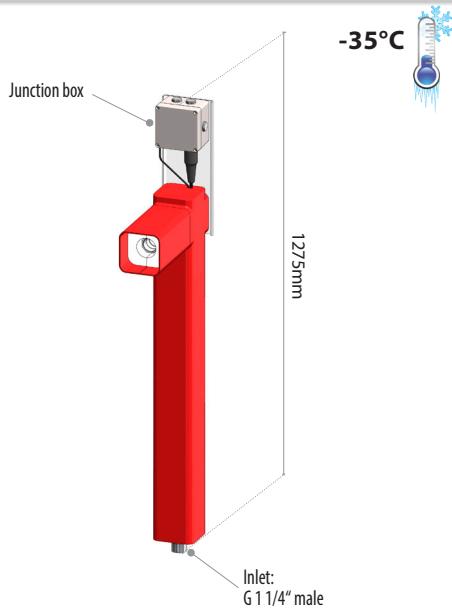


Heated underground extension pipe, trace tape heated and pre-insulated. To be used for connecting to an underground supply. Available in Galvanized or Stainless Steel.

#### TECHNICAL SPECIFICATIONS

Min. working temperature: -35°C  
Supply voltage: 230VAC, 50/60Hz

#### ORDERING INFORMATION

BROEN-LAB item no.	Material	 Ex-proof*	Min. working temperature	Mounting
H7 51 001 10 00 321	Galvanized Steel	Yes	-35°C	Underground
H7 51 001 10 00 221	Stainless Steel	Yes	-35°C	Underground

\*Ex-proof zone 1 and zone 2.

### Foot Panel for Operation of Heated Eye Wash



Weight: Approx. 23 kg

**No. 17 803 009** - Stainless Steel version

**No. 17 833 009** - Galvanized Steel version

Foot panel or treadle for operation of **heated body shower** is also available - contact your BROEN-LAB sales representative for further information.

### Foot Treadle for Operation of Heated Eye Wash



**No. 17 804 009** - Stainless Steel version

For heated eye washes operational down to -20°C

**No. 17 805 009** - Stainless Steel version

For heated eye washes operational down to -35°C

For more information and technical datasheets please contact your BROEN-LAB sales representative.



### Emergency Showers Are Important

Thousands of industrial accidents happen every year in spite of precautionary safety measures, accidents that are caused by human exposure to heat (fire, boiling water and oil, steam, etc.), to corrosive agents (fluid acids and caustics) and/or to influence by dust and corrosive aerosols (metal, dirt, fluids).

These accidents cause the injured person acute pain, disablement and may even lead to loss of life.

However, experience shows that the consequences and human suffering following such accidents can be considerably limited if water is applied directly to the affected areas instantly after the exposure.

It is a fact of life that the risk lurks in many places: in the laboratory, in the kitchen, in the warehouse, in the chemical plant, in the workshop or generally speaking in any place where heat, corrosion or dust can cause serious human injury.

For these reasons we strongly recommend the installation of adequate emergency shower systems as a means of limiting these risks.

### Location and Operation of Emergency Showers

As a basic rule the shower must be located as close to the workplace as possible so that the general first aid instructions of flushing the affected area immediately after the accident has occurred can be followed.

Please note that if strong acids or caustics are involved immediate action is demanded.

When installing the emergency shower system please also be aware that it must be easy to locate and access at all times for potential users.

The BROEN-LAB showers can be recognized by their red body shower pull rods and red eye shower push handles. In addition the internationally adapted graphic symbols for the clear identification of emergency shower installations are always included in BROEN-LAB solutions.

Further information and guidelines can be found in the EN 15154 part 1 & 2 - European standard for Emergency safety showers

### Operation

All BROEN-LAB emergency shower systems are designed and engineered to secure easy operation. This allows for "one-action" activation of the shower system ensuring immediate and continued flushing without any need for further action.

This easy operation also allows the injured person to use both hands to keep his eyes open for effective flushing, to undress under the shower or just to help exposing the affected area for better and more expedient treatment.

On the eye shower the nozzles are protected by rubber cups to keep the injured person from hurting himself unintentionally. The nozzles are also protected against dust and dirt by the caps which are removed automatically by the initial jet of water.

### Water Pressure and Temperature

Emergency showers shall be connected to a cold water system with at least drinking water quality which ensures clean and fresh water. In addition we recommend the water temperature to be kept within the range of 15° to 25° Celsius in order to provide the best immediate effect, reasonable human comfort and at the same time avoid hypothermia.

All eye showers are delivered with FLOWFIX water regulating devices ensuring constant adequate water flow regardless of high water pressure.

The illustrations will inform you of the capacities of the showers.

### Threading

In general the BROEN-LAB Emergency Showers are equipped with threadings as follows:

E.g. G 1/2": In accordance to ISO 228/1.

E.g. NPT 1/2": In accordance to ANSI-B1.20.1.

### Maintenance and Information

When installing emergency shower systems, the relevant safety advisor must provide instructions on operation and maintenance!








As basic information the potential user must be advised of:

- The location of the emergency shower system
- How to operate it and
- How to behave in case of an emergency

BROEN-LAB recommends users to check installed emergency shower systems once a week.

### IN CASE OF AN EMERGENCY

In the event of accidents caused by human exposure to heat, corrosive agents and/or to influence of dust and corrosive aerosols:

-  Start flushing immediately – the seconds count.
-  Remove contaminated clothing which is not stuck to the skin.
-  Continue flushing.
-  In case of an eye injury keep the eyelid open when continuing flushing.
-  Call for an ambulance. Ask the doctor for advice.
-  Continue flushing in the ambulance.
-  Keep the flushing in the emergency room and during transport to the special ward.

In case of blistering do inform the doctor and the hospital about the substance in question and about possible countermeasures.

RAPID ACTION CAN BE OF VITAL IMPORTANCE.

### PERSONAL WASH UNIT

The first seconds following an eye injury are often critical to keeping eye injury to a minimum. A personal wash unit may be kept in the immediate vicinity of employees working in a potentially hazardous area.

The main purpose of these units is to supply immediate flushing. With this accomplished, the injured individual should then proceed to a plumbed or selfcontained eyewash and flush the eyes for the required 15-minute period.

REFERENCE: A1 PERSONAL WASH UNIT, ANSI Z358.1.

# INTERNATIONAL STANDARDS

## NORM EN 15154-1&2


(European Committee for Standardization)

### MAIN REQUIREMENTS


- The water supplied by the body shower shall be of constant flow rate in accordance with national regulations at a flow pressure to be specified by the manufacturer. The flow pressure shall be specified and measured where the shower is connected to the water system. The body shower shall be capable of delivering this supply for a minimum of 15 minutes. Note: When no national or local regulations apply, a constant flow rate of min. 60 l/min is suitable.
- It shall only be possible to make adjustments with a tool to the direction of spray or the water distribution of a shower head. The shower head shall be self-draining between the valve and the outlet. The shower head shall be removable for maintenance but only by use of a tool.
- At a distance 700 mm below the shower head, 50% (+/-10%) of the volume of water delivered shall fall in a circle with a radius of 200mm; the water level in the individual compartments in this circle shall not deviate by more than 30% from the mean value. At this measuring level, the area reached by minimum 95% of the water shall be limited to a circle with a radius of 400 mm. The velocity of the water spray shall be low enough to be non-injurious to the user.
- The shower head shall be designed to be installed so that its lower edge is (2200 +/- 100) mm above the level on which the user stands.
- The valve shall not close automatically once it has been opened. The direction of operating the valve actuator shall be clearly visible and unmistakable.

Flow rate: Min. 60 l/min


### OPEN SPACE

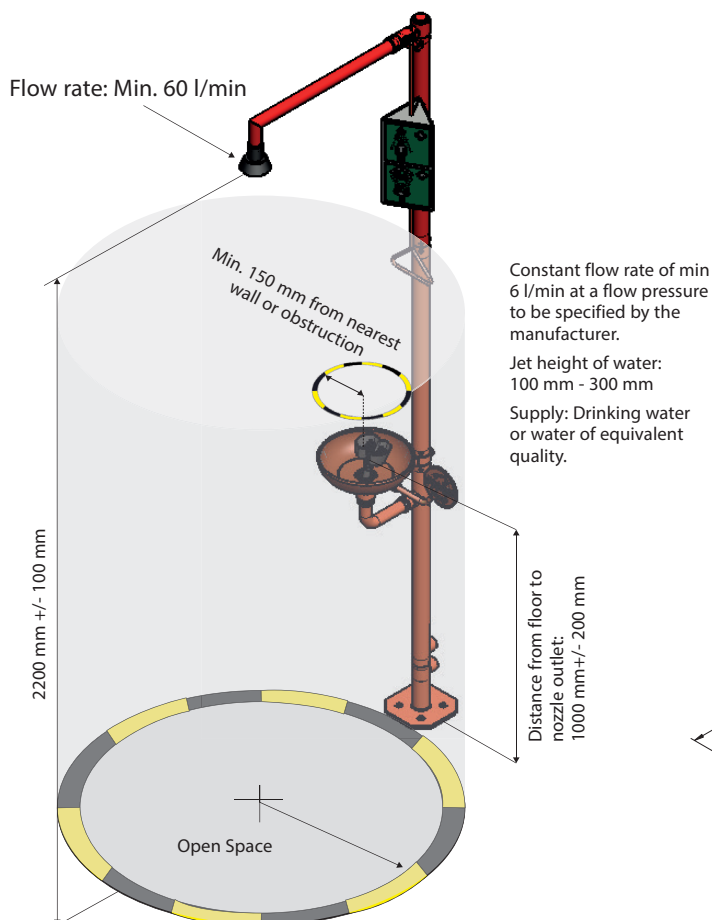

### WATER PATTERN

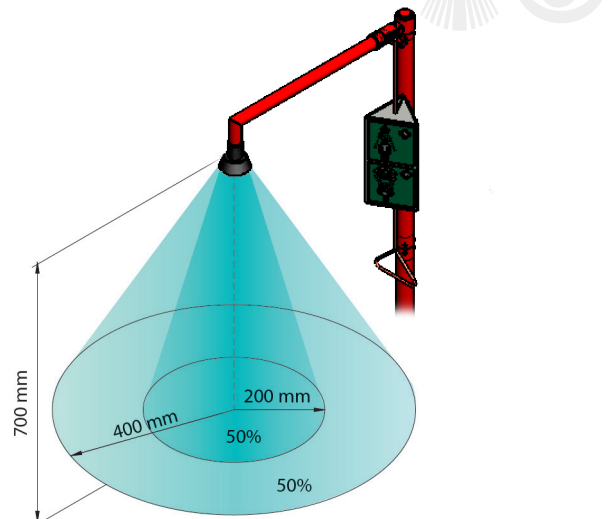

### TYPE TESTING OF WATER DISTRIBUTION

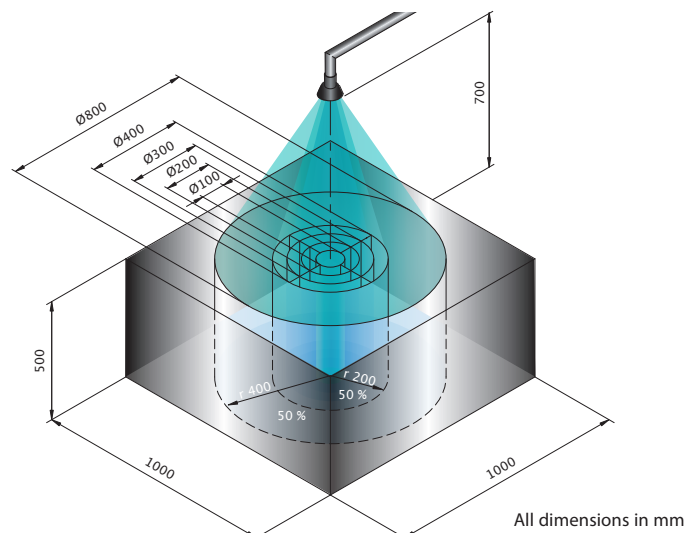



Please refer to [www.BROEN-LAB.com](http://www.BROEN-LAB.com) for a 3rd party certificate for EN 15154

(American National Standard Institute)

- Emergency eyewash and showers must be located in such a way that they can be reached by the user within maximum 10 seconds and available for immediate use (ref.: 4.5.2 + Safety considerations - B5).
- Eye/face showers must be able to deliver at least 11.4 l. water per minute continuously for at least 15 minutes (ref.: 6.1.6).
- The handle of pull rod must be located at a maximum height of 173.3 cm above the floor (ref.: 4.2).
- The valve must go from "off" to "on" position within 1 second or less and be able to remain open without further action. (ref.: 4.2).
- The body shower must be capable of delivering flushing fluid at a minimum of 75.7 liters per minute continuously for a minimum of 15 minutes (ref.: 4.1.4).

- The diameter of the water spray pattern must be minimum 50.8 cm measured at 152.4 cm above the floor (ref.: 4.1.4).
- The face of the showerhead shall be located within a range of 208.3 cm to 243.8 cm from the floor (ref.: 4.4.1).
- The eye/face showers flushing fluid flow pattern must be located within a range of 83.8 cm to 134.6 cm above the floor (ref.: 5.4.4).
- The shower must be activated at least once a week (ref.: 4.6.2).
- Tepid water is required in all showers (ref.: 4.5.6 + Safety considerations - B5).

Reference: ANSI Z 358.1-2014

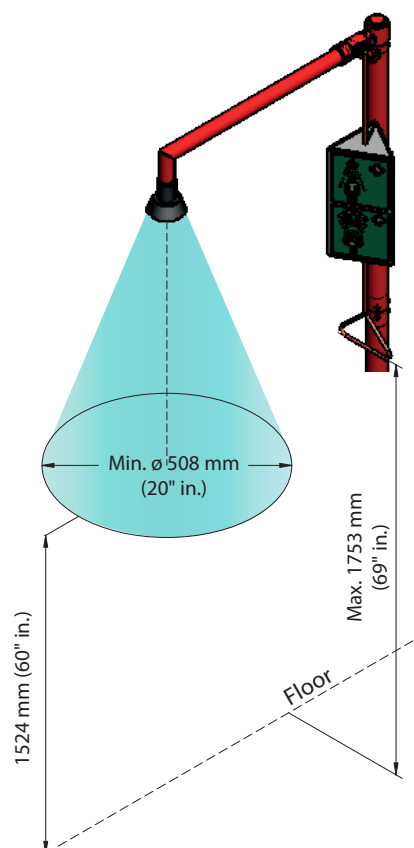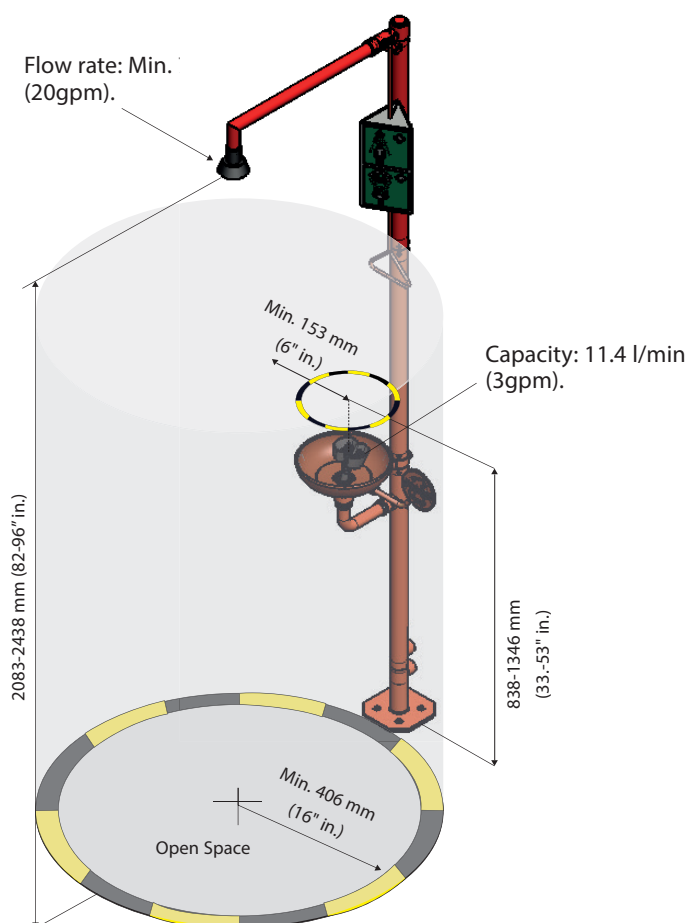

## EMERGENCY SHOWER SYSTEMS

*An eye shower is the only allowable exception from the rule that no obstructions are allowed in the "open space".*

## SPECIFICATION OF MATERIALS

Shower outlet	Brass/PP or PPO
Stand pipe	Stainless steel (AISI 304) or brass
Rubber cups	EPDM
Dust caps	POM
Bowl	Stainless steel (AISI 304)

Available in chemical resistant "BROEN-LAB Polycoat" (polyester powder coating) or "Stanniloy" plating (a highly anti-corrosive treatment based on a tin/nickel alloy).

## TEKNISKE SPECIFIKATIONER

**Min. working pressure (Dynamic/Flow pressure):** See individual products.

**Recommended max. working pressure:**

kPa	kN/m <sup>2</sup>	bar	psi
500	500	5	72.5


**Max. working pressure (Static):**

kPa	kN/m <sup>2</sup>	bar	psi
1000	1000	10	145

**Recommended temperature:** 25°C.

Alternatively, install TMV.

## FLOWFIX

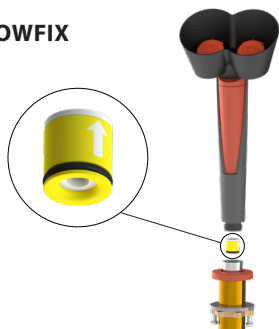

FLOWFIX is a feature that is built into most of BROEN-LAB's eye baths and hand/eye showers. Flowfix is a water regulating device, that controls the water flow at a preset rate regardless of water pressure.

## SYMBOLS:



Eye shower



Body shower




For disabled people as well

## SELF-DRAINING SHOWER HEAD


BROEN-LAB has developed a new self-draining emergency shower head in compliance with both EN15154-1 and ANSI Z358.1.

BROEN-LAB has 3rd party certificate for EN 15154 and ANSI Z358.1-2014. The self-draining shower head requires min. 1 bar (14.5 psi) and max. 10 bar (145 psi) dynamic water pressure for optimal function.

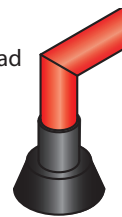

## FEATURES & BENEFITS

- ✓ Self-draining, anti lime scale - prevents bacteria build-up and reduces maintenance
- ✓ Complies to both EN & ANSI Standards
- ✓ Chemical resistant PPO plastics
- ✓ Soft water spread
- ✓ Available for 1/2", 3/4" and 1" pipes

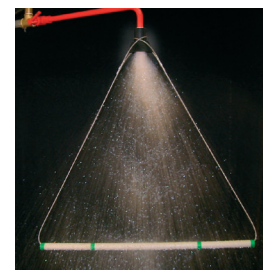

## CALCULATING THE FLOWRATE AND DYNAMIC PRESSURE DROP

Any calculation can be done using the following formulas:

$$K_v = \frac{q_v [\text{m}^3/\text{h}]}{\sqrt{\Delta P [\text{bar}]}} \quad \text{or} \quad q_v = K_v \times \sqrt{\Delta P} \quad \text{or} \quad \Delta P = \left( \frac{q_v}{K_v} \right)^2$$

## Where

$K_v$  = flowrate [m<sup>3</sup>/h] or [l/sec] at 1 bar pressure drop

$q_v$  = flowrate [m<sup>3</sup>/h] or [l/sec]

$\Delta P$  = pressure drop [bar]

1) Please note that eye showers with built-in FLOWFIX have a fixed flowrate.

## Examples:

What pressure drop is needed in order to achieve 75 l/min = 1.25 l/sec through a body shower with 1/2" inlet?

$$\Delta P = \left( \frac{q_v}{K_v} \right)^2 = \left( \frac{1.25 [\text{l/sec}]}{1.08 [\text{l/sec}]} \right)^2 = 1.34 [\text{bar}]$$

What flowrate can be expected through a free standing body shower with top inlet at a pressure drop of 1.5 bar?

$$q_v = K_v \times \sqrt{\Delta P} = 1.10 [\text{l/sec}] \times \sqrt{1.5 [\text{bar}]} = 1.35 \text{ l/sec} = 80.8 \text{ l/min.}$$

## K<sub>v</sub> VALUES

Type of shower	Eye shower without FLOWFIX <sup>1)</sup>	Body shower 1/2" pipe or inlet	Body shower 3/4" pipe or inlet	Body shower 1" pipe or inlet	Free standing-body shower, top inlet	Free standing-body shower, bottom inlet	Thermostatic mixing valve 17 140.009	Thermostatic mixing valve 17 512 009/ 170513 009	
K <sub>v</sub> -value	m <sup>3</sup> /h	0.75	3.89	4.3	5.31	3.97	2.99	1.25	4,7
	l/sec	0.21	1.08	1.19	1.48	1.10	0.83	0.35	1.33



This image shows a blank sheet of white paper designed for note-taking. At the top, there is a solid red horizontal bar. Below this bar, the word "NOTES" is printed in a bold, black, sans-serif font. The remainder of the page is filled with thin, light gray horizontal lines spaced evenly apart, providing a guide for writing.

NOTES

## This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There is no text or other markings on the paper.

BROEN-LAB develops, manufactures and sells Laboratory Fittings, Emergency Shower Systems and Eye Wash Systems into a broad spectrum of laboratories and industrial locations; sectors include Pharmaceutical, Food&Beverage and Academia. Our expertise and product quality ensure optimal solutions compliant to all relevant international norms setting new standards in risk mitigation in modern research and development facilities.

BROEN-LAB offer solutions that ensure the functionality, compliance, hygiene, durability and safety are of the highest priority. For over 50 years our products have been integrated into a wide variety of workplaces, laboratories, hospitals and industrial locations, with features that are assessed and designed to be flexible, durable and compliant offering a broad range of options tailored to each location. This inherent design integrity provides confidence in the solution adopted wherever in the world these are applied.

BROEN-LAB is a collaboration partner in the design and layout of your next laboratory and it's integrated safety

## HEADQUATER

### **BROEN-LAB A/S**

Drejervaenget 2 · DK-5610 Assens  
Tel. +45 6376 6376 · Fax +45 6471 2476  
E-mail: lab@broen-lab.com



### **BROEN-LAB A/S ISO 9001 certification**

In September 1991 BROEN-LAB A/S was certified according to ISO 9001 as one of the first Danish companies. The certification was carried out by Bureau Veritas, Denmark, for our Danish site in Assens. The quality management system of BROEN-LAB A/S now complies with detailed specifications laid down by the internationally acknowledged EN ISO 9001:2015. This certification will further contribute to reputation for quality and reliability of BROEN-LAB A/S.

## INTERNATIONAL OFFICES

### **BROEN-LAB Ltd**

Tel. +44 121 522 4515  
E-mail: lab@broen-lab.com

### **BROEN-LAB Singapore Pte. Ltd**

25 International Business Park ·  
#04-60A German Centre · Singapore 609916  
Tel. +65 6298 0662 · FAX +65 6298 0468  
E-mail: lab@broen-lab.com

### **BROEN-LAB GmbH**

Tel. +45 6376 6376 · Fax +45 6471 2476  
E-mail: lab@broen-lab.com

### **BROEN-LAB Sverige AB**

Tel. +45 6376 6376 · Fax +45 6471 2476  
E-mail: lab@broen-lab.com

broen-lab.com

**BROEN-LAB**